
Speed Read Digital Thermometer

Model 78-359

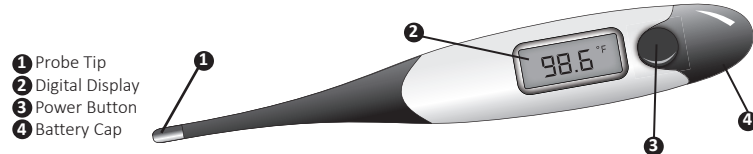THANK YOU FOR PURCHASING A WELBY SPEED READ DIGITAL THERMOMETER. PLEASE READ THESE INSTRUCTIONS CAREFULLY TO ENSURE PROPER AND SAFE OPERATION OF YOUR THERMOMETER.

PLEASE READ AND SAVE THESE INSTRUCTIONS

- This product is not intended to diagnose or treat any health problem or disease and should not be used as a substitute for the consultation and advice of a physician or other medical professional.
- Thermometers should only be used under the supervision of an adult.
- This thermometer should be held during use; the patient should be seated or lying down during readings.
- Before use, it is recommended that the patient avoid eating, drinking, smoking or exercising for at least twenty minutes prior to taking a measurement. Never walk, run or talk during a measurement.

NORMAL BODY TEMPERATURE AND FEVER

Body temperature can vary from 97°F (36.1°C) and 99°F (37.2°C) and still be considered normal. The medically accepted 'normal' body temperature is 98.6°F (37°C) Body temperature is commonly lower in the morning than at any point during the rest of an individual's waking hours.

Temperature results will vary based on the body location point of the reading. Oral temperature readings between 97°F (36.1°C) and 100°F (37.7°C) are considered normal. A rectal temperature reading is generally 1°F (0.5°C) higher; while an underarm, or axillary, temperature will be 1°F (0.5°C) lower.

A fever is defined as a body temperature that is elevated above the normal for that person. It is important to determine what is normal for an individual before determining if a fever is present. Tracking an individual's temperature on a consistent basis, at the same body site and at the same time every day will help determine an individual's normal temperature.

FAHRENHEIT CELSIUS SWITCHABLE FEATURE

This unit takes measurements in both Fahrenheit and Celsius. While the unit is OFF, press and hold the POWER button. After the brief self-test, the unit will show the last reading in either °F or °C. Continue to hold the POWER button, and display will change the mode and display the new temperature reading mode. Release the POWER button and begin your reading.

MEMORY RECALL

This unit includes a recall of the last measurement taken. When the unit is OFF, press the POWER button to turn the unit on, and release. After the self-test, the unit will briefly show the last reading with an 'M' on the display.

OPERATING INSTRUCTIONS

Press and release the POWER button to activate the thermometer. The unit will beep and the display will show all of the digital display symbols as it performs a self-test before showing the last reading in memory.


A Lo with a flashing °F will appear in the display, followed by a brief flash of 98.6°F (37°C); the display will return to the Lo screen, indicating that the thermometer is ready for measurement.

ORAL USE

Approximate measurement time: 10-20 seconds


- Apply a probe cover if desired.
- Place the probe under the tongue as near as possible to a heat pocket as noted in the diagram. Readings will be progressively lower the further the probe is from the indicated heat pockets.
- Hold the thermometer in place during the reading; do not bite down on the thermometer. The mouth must remain closed to ensure accurate results.
- When the reading is complete and the peak temperature has been reached, the °F will cease flashing. If the reading is below 100°F (37.7°C), the unit will sound ten long beeps; if the reading is at or above 100°F (37.7°C), the unit will sound ten sets of three short beeps.
- To turn the unit off, simply press the POWER button; the unit will automatically shut-off after approximately ten minutes.

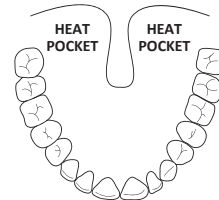UNDERARM (AXILLARY) USE

Approximate measurement time: 20-30 seconds

- Wipe the underarm with a dry towel before application. Apply a probe cover if desired.
- Place the probe tip under the arm so that the tip is touching skin; hold the thermometer parallel, slightly angled upwards as indicated in the diagram. Position the arm across the chest to ensure the tip is completely covered.
- When the reading is complete and the peak temperature has been reached, the °F will cease flashing. If the reading is below 100°F (37.7°C), the unit will sound ten long beeps; if the reading is at or above 100°F (37.7°C), the unit will sound ten sets of three short beeps.
- To turn the unit off, simply press the POWER button. The unit will automatically shut-off after approximately ten minutes.

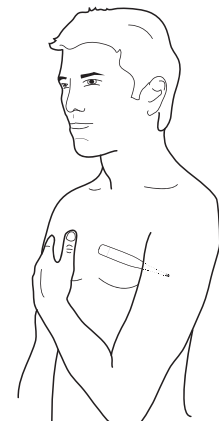Axillary temperature readings are generally 1°F lower than oral readings.

RECTAL USE

Approximate measurement time: 9 seconds

- Apply a probe cover if using; you may use water-soluble lubricant for easier insertion if desired.
- Lay the patient on their side or stomach with their legs hanging down.
- Gently insert the tip of the probe NO MORE THAN 1 INCH into the rectum. Stop immediately if you detect resistance. Always hold the thermometer in place during measurement.
- When the reading is complete and the peak temperature has been reached, the °F will cease flashing. If the reading is below 100°F (37.7°C), the unit will sound ten long beeps; if the reading is at or above 100°F (37.7°C), the unit will sound ten sets of three short beeps.
- To turn the unit off, simply press the POWER button. The unit will automatically shut-off after approximately ten minutes.

Rectal temperature readings are generally 1°F higher than oral readings.

PROBE COVER USE

Disposable, hygienic probe covers are recommended for use with stick thermometers to prevent contamination. Probe covers are intended for one-time use only; a new probe cover should be used with each application. Probe covers can commonly be found at most local pharmacies and drug store retailers.

CARE AND CLEANING INSTRUCTIONS


Store the unit in the provided protective case when not in use. Do not store the thermometer where it will be exposed to direct sunlight, dust or humidity. Avoid extreme temperatures. Do not attempt to disassemble the unit. If ERR displays and no measurement can be obtained, the unit must be replaced.

Clean the thermometer by wiping with rubbing alcohol (70% isopropyl alcohol) or soap and water. Do not submerge in any liquid or autoclave. Wipe the thermometer with a dry cloth and allow to air dry.

REPLACING THE BATTERY

The low-battery indicator symbol will appear on the lower right of your display when the battery needs to be changed.

- Remove the battery cap by pulling gently.
- Looking at the unit from the back side (ratings label face up) carefully slide the battery tray out so the entire battery is exposed.
- Remove the battery with a non-metal tool. Do not attempt to remove the thin plastic battery cover bar.
- Replace the battery with a LR41 alkaline or SR41 silver oxide type battery with the positive (+) side facing up.
- Carefully replace the battery cap.
- Do not dispose of used batteries in the trash. Consult your local recycling center for the proper way to dispose of batteries in your area. Keep used batteries away from children, pets and excessive heat.

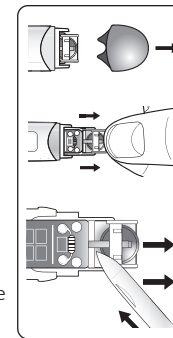SPECIFICATIONS

BATTERY:	One 1.5V button size battery LR41 or equivalent
BATTERY LIFE:	Approximately 24 months
AUTOMATIC SHUT-OFF:	Approximately after 10 minutes of non-use
ACCURACY:	±0.2°F, 98.0°F - 102.0°F (±0.1°C, 37.0°C - 39.0°C) ±0.3°F, 96.4°F - 98.0°F (±0.2°C, 35.8°C - 37.0°C) ±0.3°F, 102.0°F - 106.0°F (±0.2°C, 39.0°C - 41.0°C)
MEASUREMENT RANGE:	90.0°F - 109.9°F (32° - 42.99°C)
OPERATING TEMPERATURE:	50.0°F - 104.0°F (10°C - 40°C)
STORAGE TEMPERATURE:	14.0°F - 140.0°F (-10°C - 60°C)
DISPLAY:	Liquid crystal display
DISPLAY RESOLUTION:	4 digits (1/10th of a degree increments)
DIMENSIONS:	5-1/2" X 7/8" X 3/8"
WEIGHT:	0.38 ounces

Specifications are subject to change without notice.

TWO YEAR LIMITED WARRANTY

Congratulations on your purchase of a Digital Thermometer. Your Digital Thermometer is covered by the following limited warranty commencing upon the date of purchase, and subject to the following terms and conditions:

The Warrantor warrants that its Digital Thermometer will be free from defects in materials and workmanship under normal consumer usage for a period of one year for the original purchaser of the Digital Thermometer.

Periodic maintenance, repair and replacement of parts due to normal wear and tear are excluded from coverage. Defects or damage that result from: (a) improper operation, storage, misuse or abuse, accident or neglect, such as physical damage (cracks, scratches, etc.) to the surface of the product resulting from misuse; (b) contact with liquid, water, rain, extreme humidity or heavy perspiration, sand, dirt or the like, extreme heat, or food; (c) use of the Digital Thermometer for commercial purposes or subjecting the Digital Thermometer to abnormal usage or conditions; or (d) other acts which are not the fault of the Warrantor, are excluded from coverage. This warranty does not cover batteries or other power sources that may be provided with, or used with the Digital Thermometer.

If the Digital Thermometer fails to conform to this limited warranty, return the Digital Thermometer postage prepaid to: **Attn: Repair Department, 1175 Lakeside Drive, Gurnee, IL, 60031**. When returning a product, please also include: (i) a copy of your receipt, bill of sale or other comparable proof of purchase; (ii) a written description of the problem; and (iii) your name, address and telephone number. Carefully package the product to avoid any damage that may occur while in transit; shipping insurance with returned receipt is recommended. At our option, Veridian Healthcare will repair or replace the unit found to be defective in materials or workmanship under normal consumer usage. The purchaser will be notified of any additional repairs required prior to completing the repair, and will be responsible for parts charges, if any, and repair charges not covered by this limited warranty.

EXCEPT AS PROVIDED FOR IN THIS LIMITED WARRANTY, ALL EXPRESS AND IMPLIED WARRANTIES AND CONDITIONS ARE DISCLAIMED, INCLUDING WITHOUT LIMITATION THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. THE REPAIR OR REPLACEMENT AS PROVIDED UNDER THIS LIMITED WARRANTY IS THE EXCLUSIVE REMEDY OF THE CONSUMER, AND IS PROVIDED IN LIEU OF ALL OTHER WARRANTIES, EXPRESS OR IMPLIED. IN NO EVENT SHALL VERIDIAN HEALTHCARE BE LIABLE, WHETHER IN CONTRACT OR TORT (INCLUDING NEGLIGENCE) FOR DAMAGES IN EXCESS OF THE PURCHASE PRICE OF THE PRODUCT, OR FOR ANY INDIRECT, INCIDENTAL, SPECIAL OR CONSEQUENTIAL DAMAGES OF ANY KIND, OR FOR DAMAGES TO, OR LOSS OF, OTHER PROPERTY OR EQUIPMENT OR PERSONAL INJURIES TO THE FULL EXTENT THESE DAMAGES MAY BE DISCLAIMED BY LAW.

Some states and jurisdictions do not allow the limitation or exclusion of incidental or consequential damages, or limitation on the length of an implied warranty, so the above limitations or exclusions may not apply to you. This warranty gives you specific legal rights, and you may also have other rights that vary from state to state or from one jurisdiction to another.

AFTER SALES SUPPORT • SERVICIO POSVENTA

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ALDI Inc.