User Manual

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dentiguard

RECHARGEABLE SONIC TOOTHBRUSH

CEPILLO DENTAL SÓNICO RECARGABLE

English.....3


Español.....9

Contents

Parts List.....	3
General Information.....	4
Safety Instructions	5
User Instructions.....	6
Cleaning and Maintenance.....	7
Other Useful Information.....	8
Technical Specifications	8
Special Instructions	8
Storing	8

MAIN PARTS

1. Brush head
2. Handle
3. Charging base with storage
4. Additional brush heads
5. LED charging indicator light
6. LED mode indicator lights
7. Power button – ON/OFF
8. Brushing modes button – Whitening / Standard / Sensitive

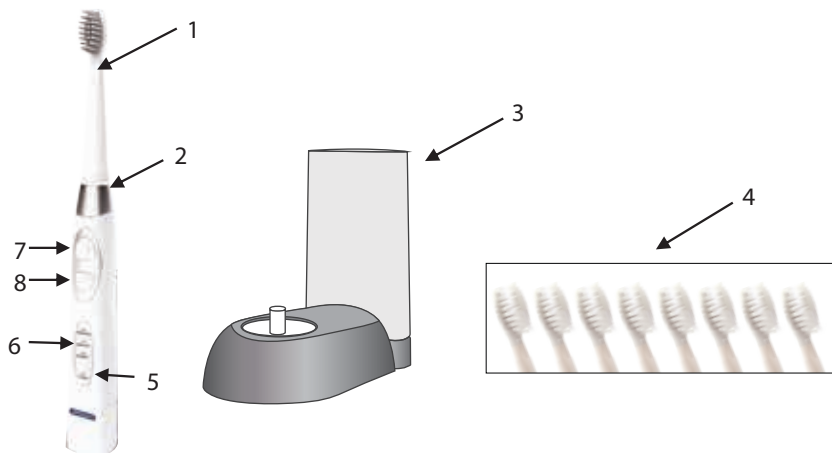**CONTENTS OF PACKAGING**

1. Charging base
2. Handle
3. 8 brush heads
4. Instruction Manual

GENERAL INFORMATION AND SAFETY INSTRUCTIONS

When using an electrical appliance, basic safety precautions should be observed, including the following.

PLEASE READ AND SAVE ALL INSTRUCTIONS BEFORE USING THE SONIC TOOTHBRUSH.

1. Always unplug the charging base after the recharging cycle is complete or if the unit will be left unattended for long periods of time.
2. DO NOT use while bathing.
3. DO NOT allow the charging base to become wet or damp.
4. DO NOT place or store product where it can fall or be pulled into a tub or sink.
5. DO NOT reach for the charging base if it has fallen into water or any other liquid. Unplug immediately.
6. Close supervision is necessary when this product is used by, on, or near children or invalids.
7. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
8. Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or if it is dropped in water. Call the customer service number provided.
9. Keep the cord away from heated surfaces.
10. DO NOT use this product if you feel drowsy.
11. DO NOT use this product outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
12. NO object(s) should ever be inserted or dropped into any opening on this product.
13. A product should never be left unattended when plugged in.
14. Never block the air openings of the product or place it on a soft surface, such as a bed or couch, where the air openings may be blocked. Keep the air openings free of lint, hair, and the like.

IMPORTANT MEDICAL SAFEGUARDS

1. Contact your dentist before using the Sonic Toothbrush if you have had oral or periodontal surgery.
2. Contact your dentist if excessive bleeding occurs or if bleeding continues after the first week of using the Sonic Toothbrush.
3. If you have questions concerning usage with a pacemaker, Please consult your physician prior to use.
4. This unit is designed for cleaning teeth and gums only. DO NOT use it for any other purposes.
5. The Sonic Toothbrush can be used with any leading brand of toothpaste.

Warning! Read and understand all instructions. Failure to follow all instructions listed below may result in electric shock, fire or serious personal injury. The warnings, cautions, and instructions discussed in this instruction manual cannot cover all possible conditions and situations that may occur. It must be understood by the operator that common sense and caution are a factor which cannot be built into this product, but must be supplied by the operator.

CAUTION: A short power supply cord is provided to reduce the risk of personal injury resulting from becoming entangled in or tripping over a longer cord.

ELECTRIC POWER

If electric circuit is overloaded with other appliances, your appliances may not operate properly. The appliances should be operated on a separate electrical circuit from other operating appliances.

SAFETY WARNING

Although your appliance is easy to operate, for your safety, the warnings below must be followed:

1. Do not place the appliance near a heat source.
2. Use the appliance on a stable surface.
3. Do not unplug by pulling on the supply cord.
4. The appliance must be unplugged:
 - A. Before any cleaning or maintenance.
 - B. If it appears to be faulty.

3-PRONG POLARIZED PLUG

NOTE:

This Sonic Toothbrush has a polarized plug (one blade is wider than the other).

As a safety feature to reduce the risk of electrical shock, this plug is intended to fit in a polarised outlet only one way.

If the plug does not fit fully in the outlet, reverse the plug.

If the plug still does not fit, contact a qualified electrician.

Do not attempt to defeat this safety feature.

THIS PRODUCT IS INTENDED FOR HOUSEHOLD USE ONLY

SAVE THESE INSTRUCTIONS

USER INSTRUCTIONS

1. Place the Sonic Toothbrush base on a level and dry surface. Plug the power cord into a polarized 100 - 240V AC wall outlet.
2. Place the handle into the base. The LED charging indicator light will light up on the handle to indicate the handle is charging. Please note the LED indicator light will stay on so long as the handle is sitting in the charging base and it is plugged in.
3. It will take approximately 24 hours for a full charge the first time. You may unplug the charger at that time if you prefer. The charger may stay plugged in but we recommend that you unplug it if you will be away for an extended period of time, i.e. vacation.
4. You should fully discharge the handle at least every 6 months in order to maintain the maximum capacity of the rechargeable batteries. You can achieve this by unplugging the base and letting the handle run down with normal use.
5. Insert the brush head onto the shaft of the handle. There will be a small space (1mm) between the brush head and the handle. This allows for a free motion of the brush head during brushing.
6. Apply toothpaste to the brush head before switching it on.
7. Select the desired mode, Whitening, Standard, Sensitive by repeatedly pushing the power button. The mode indicator light will light up with each push and will switch and then OFF.
8. Move the brush head gently along your teeth and gum line and allow the brush heads to stay on the surface of each tooth for a few seconds with a slight circular motion for the best cleaning and gum massaging effect.
9. This unit comes with a built in timer that will pulse 3 times at the end of 30 seconds. This allows you to concentrate on each quadrant of your mouth for 30 seconds with a total brushing time of 2 minutes. DO NOT forget to reach all areas of your teeth including the inner and chewing surfaces and behind your back teeth.
10. After 2 minutes the handle will automatically turn off. After use, rinse the brush head under a tap and wipe the handle dry.

NOTE: When you first use your Sonic Toothbrush you may experience slight bleeding. This should stop after a few days when your gums get used to the toothbrush. If the bleeding persists, please consult your dentist or hygienist.

CLEANING AND MAINTENANCE

1. Disconnect the charging base cord from the outlet.
2. DO NOT immerse the base, power cord or plug in water or any other liquids.
3. Brush Heads: After each use, rinse the brush head for several seconds under running water. You may soak the brush head in denture cleanser every 2 to 4 weeks to keep it hygienically clean.
4. Charging Base and Handle: Periodically clean the exterior of the handle and charging base with a soft slightly damp towel or paper towel. Dry thoroughly before plugging back into the outlet and placing the handle back on the charger.
5. DO NOT wrap the power cord around the charging base.
6. This unit has no serviceable parts, therefore NO altering or servicing should be performed by any unauthorized representative.

REMOVABLE RECHARGING BASE

This unit can be removed for easier washing if you prefer.

OTHER USEFUL INFORMATION

Technical Specifications

Model Number	46365
Rated Wattage	0.9 – 1.5W
Rating Voltage	100 - 240V
Rechargeable battery	Lithium Battery

Special Instructions

The Sonic Toothbrush is equipped with a long- lasting rechargeable Lithium battery. To help protect our environment, it should be removed at the end of the unit's life for safe disposal.

THE BATTERY CANNOT BE REPLACED. To remove the battery before discarding the unit:

1. Open the battery cover at the bottom of the handle using a Phillips screwdriver to remove the screw. Remove the end cap.
2. Remove the lithium battery. You may need to gently shake the handle to dislodge it.
3. Dispose of the battery in compliance with your local regulations for environmental protection and recycling of raw materials.

Storing

Store the unit in a dry location. Do not place any heavy items on top of unit during storage as this may result in possible damage of Sonic Toothbrush.

Service Center

If you have any questions in regards to the operation of this appliance please contact our service center at:

Ashtel Studios Inc.

Tel: 1-877-274-8358

(Business hours Mon – Fri 11:00 AM – 8:00 PM EST)

Email: support@brushbuddies.com

Environmental Protection

If the appliance should no longer work at all, please make sure that it is disposed of in an environmentally friendly way.

FCC STATEMENT

NOTE: This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.
- This device complies with part 18 of the FCC Rules.

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ALDI INC. FULL WARRANTY CONDITIONS

Dear Customer,

The **ALDI warranty** a full warranty offering you the following benefits:

Warranty period: 2 years from date of purchase.

Costs: Free repair/ exchange.
No transport costs.

Hotline: (877) 247-8358

ADVICE: Please contact our service hotline by phone, e-mail or fax before sending the device. This allows us to provide support in the event of possible errors.

In order to make a claim under the warranty, please send us:

- The faulty item together with the original purchase receipt and the completed warranty card.
- The product with all components included in the packaging.

The warranty does not cover damage caused by:

- **Accident** or **unanticipated events** (e.g. lightning, water, fire).
- **Improper use** or **transport**.
- **Disregard of the safety** and **maintenance instructions**.
- Other **improper treatment** or **modification**.

After the expiry of the warranty period, you still have the possibility to have your product repaired at your own expense. If the repair or the estimate of costs is not free of charge you will be informed accordingly in advance.

This warranty does not affect your statutory rights. The period of warranty can only be extended in accordance with a legal standard. In countries where a (compulsory) warranty and/or spare part storage and/or a system for compensation are/is required by law, the statutory minimum conditions apply. In the event that a product is received for repair, neither the service company nor the seller will assume any liability for data or settings possibly stored on the product by the customer.

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