User Manual

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ViSAGE[®]

**SHIATSU FOOT MASSAGER
WITH HEAT**

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
Español.....11

Contents

Parts list	3
Contents	3
Assembly	4
General information	5
Safety instructions	7
User instructions	8
Cleaning and maintenance	9
Other useful information	10
Technical specifications	10
Storing	10

Main parts

1. Power ON/OFF button
2. Heat ON/OFF button
3. Main unit
4. Power A/C adapter
5. Power D/C plug
6. Power D/C receiver
7. Handle

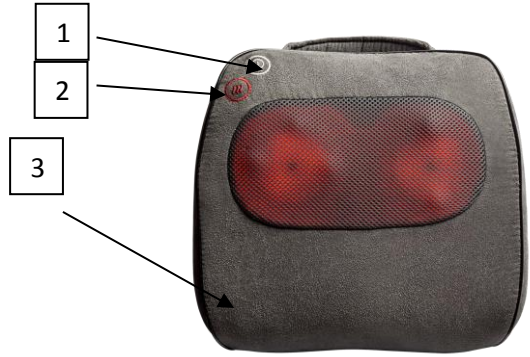**Contents of packaging**

- Shiatsu foot massager with heat
- Power adapter
- Instruction manual

Assembly

1. Unpack the Shiatsu Foot Massager from the box.
2. Remove any unwanted packaging materials such as cardboard, plastics or styrofoam.
3. Connect the power adapter plug to the power receiver on the massager.
4. Plug into outlet.

General information

When using the Shiatsu Foot Massager with Heat, basic safety precautions should always be followed. **PLEASE READ AND SAVE ALL INSTRUCTIONS BEFORE USING THE SHIATSU MASSAGER.**

- **THIS ITEM CONTAINS A HEATING ELEMENT.** Although it only provides a soothing warmth, people who are sensitive to heat should pay attention when using.
- Connect the massager to a properly grounded outlet only.
- This massager is intended for indoor, non-commercial, household use only. Do not use the massager for any other purpose. Do not use outdoors.
- Do not leave the massager unattended while it is in use. Always unplug the massager from the electrical outlet when it is not in use.
- Keep the massager out of the reach of children and pets. This massager is not intended for use for infants, children or disabled persons without close supervision.
- Use this massager only as described in this manual. Any other use not recommended by the manufacturer may cause fire, electric shock, or personal injury.
- Do not use the massager to supplement medical treatment.
- Never locate or store the massager where it may fall or be pulled into a bathtub, sink or other container with liquid in it. Keep the power supply cord and the massager away from heated areas.
- Do not reach for an appliance that has fallen into water. Unplug it immediately.
- Do not immerse in water or other liquids.
- Do not operate the massager in wet or moist conditions.
- Do not allow the power supply cord to hang (e.g., over the edge of a table or counter) or place/run the power supply cord under rugs, carpeting, or in high-traffic areas where it may be tripped over or pulled.
- Do not attempt to plug or unplug the massager while massaging your feet.
- To disconnect the massager, remove the plug from the electrical outlet. Never pull on the power supply cord when disconnecting from power outlet.
- Do not operate, plug or unplug the massager from the electrical outlet with wet hands.

- Do not operate any appliance with a damaged cord or plug, after the appliance malfunctions, or has been dropped or damaged in any manner. Do not use the massager while bathing or showering.
- Do not use attachments not recommended by the supplier, specifically any attachments not provided with the massager.
- Do not stand on the massager.
- Do not fall asleep while using the massager. Never use while sleeping. Do not use the massager before going to bed. The massager has a stimulating effect and may delay sleep.
- Never use pins or other metallic fasteners with the massager. Never drop or insert any object into any opening.
- Do not attempt to repair or disassemble the massager. There are no user serviceable parts.
- Do not operate massager where aerosol (spray) or similar products are being used.
- Do not operate the massager under a blanket or a pillow. Excessive heating can occur and cause fire, electric shock or injury to user.
- If you have concerns regarding your health, consult your doctor before using the massager. Never use the massager if you are sensitive to heat, have sensitive skin, open wounds, discolored areas, or any area of the body that is swollen, burned, inflamed or where skin eruptions or sores are present.
- Massage and heat should be pleasant and comfortable. If you experience muscle or joint pain for a prolonged period of time or if pain or discomfort should result, discontinue use and consult your doctor.
- In case of pregnancy, diabetes, or illness consult your doctor before use.
- Individuals with pacemakers should consult their doctor before using the massager.
- Keep cord away from heated surfaces.
- Do not carry the massager by its power cord or use the cord as its handle. Do not pull or twist the power cable.
- To disconnect, turn off all controls, then remove the plug from the outlet.

SAVE ALL INSTRUCTIONS

Safety instructions

WARNING: Read and understand all instructions. Failure to follow all instructions listed below may result in electric shock, fire or serious personal injury. The warnings, cautions, and instruction discussed in this instruction manual cannot cover all possible conditions and situations that may occur. It must be understood by the operator that common sense and caution are a factor which cannot be built into this product, but must be supplied by the operator.

CAUTION: A short power supply cord is provided to reduce the risk of personal injury resulting from becoming entangled in or tripping over a longer cord. Extension cords are available from local hardware stores and may be used if care is exercised in their use. If an extension cord is required, special care and caution is necessary. Also the cord must be: (1) marked with an electrical rating of 125V, and at least 13 A., 1625 W., and (2) the cord must be arranged so that it will not drape over the countertop or tabletop where it can be pulled on by children or tripped over accidentally.

Electric power

If electric circuit is overloaded with other appliances, your massager may not operate properly. The massager should be operated on a separate electrical circuit from other operating appliances.

Safety warning

Although your massager is easy to operate, for your safety, the warnings below must be followed:

1. Do not place the massager near a heat source.
2. Never leave the massager unattended while connected to the electrical outlet.
3. Do not unplug the massager by pulling on the power supply cord.
4. The massager must be unplugged:
 - A. Before any cleaning or maintenance.
 - B. After use.
 - C. If it appears to be faulty.

**THIS PRODUCT IS INTENDED FOR HOUSEHOLD USE ONLY.
SAVE THESE INSTRUCTIONS!**


User Instructions

How to use the Shiatsu Foot Massager with Heat

1. Connect the A/C power adapter plug into the D/C connector of the massager.

2. Plug the power supply cord into an electrical outlet. Be sure plug fits tightly in the outlet.
3. After washing feet and drying thoroughly, place the massager on the floor near a chair or sofa, sit down in a comfortable position and place your feet on the massager.
4. The Power button and Heat button are located in the upper left hand corner.

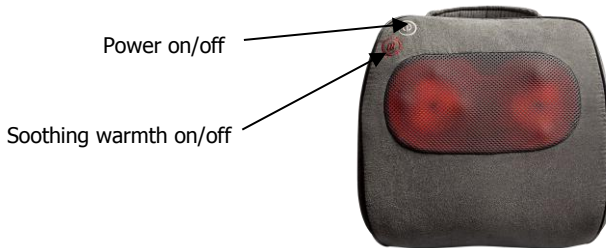5. Press the white power button to turn on. The massage heads will rotate outward.

6. Press the power button again, wait 2 seconds and the massage heads will rotate inward.

7. If you would like to add soothing warmth to the massage press the Red heat button to turn it ON. Press the Red heat button again to turn OFF the warmth.
8. The shiatsu massage can be enjoyed for as long as 15 minutes before it will automatically turn itself off. If you do not want to use the massager for 15 minutes, simply switch off the power button.

Cleaning and maintenance

CAUTION:

Always be sure that the Shiatsu Foot Massager is switched OFF, unplugged and completely cooled before cleaning or attempting to store the item. NEVER IMMERSE in water or any other liquid.

1. This massager can only be cleaned by using a soft towel or cloth.
2. NEVER use brushes, aggressive cleaning agents, gasoline, paint thinners or alcohol.

DO NOT IMMERSE THE MASSAGER IN WATER AT ANY TIME

DO NOT USE ANY ABRASIVE MATERIALS TO CLEAN SURFACE

Other useful information

Technical specifications

Model Number	TM-211
Rating Voltage	120V, 60Hz
Nominal Power	18W

Storing

- Always unplug the massager before storing.
- Store the massager in a dry location in the original box if possible.
- Do not place any heavy items on top of the massager during storage as this may result in possible damage.

Service center

If you have any questions in regards to the operation of this massager or are in need of a spare part, please contact our service center at:

Tel: 1-888-367-7373

Business Hours: Mon-Fri 10:00am – 5:00pm EST

Email: help@myproduct.care

Environmental protection

If the massager should no longer work at all, please make sure that it is disposed of in an environmentally friendly way.

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 _____ Email _____

Date of purchase* _____

*We recommend you keep the receipt with this warranty card

Location of purchase _____

Description of malfunction:

Return your completed warranty card to:

Wachsmuth & Krogmann, Inc.

1015 Hawthorn Drive

Itasca, IL

60143

USA

help@myproduct.care

AFTER SALES SUPPORT

49581

 888 367 7373

 help@myproduct.care

MODEL:

TM-211

05/2017

Hotline:

888-367-7373

Operating hours:

Monday – Friday

10:00am – 5:00pm EST

2

YEAR WARRANTY
AÑOS DE GARANTÍA

ALDI INC. FULL WARRANTY CONDITIONS

Dear Customer:

The **ALDI warranty** is a fully warranty offering you the following benefits:

Warranty period: **2 years** from date of purchase.
 6 months for movable parts and consumables under normal and proper conditions of use (e.g. rechargeable batteries).

Costs: Free repair/exchange or refund.
 No transport costs.

ADVICE: Please contact our service hotline by phone, e-mail or fax before sending in the device. This allows us to provide support in the event of possible operator errors.

In order to make a claim under the warranty, please send us:

The faulty item with all the component parts, the original receipt and the warranty card properly completed.

The warranty does not cover damage caused by:

Accident or unanticipated events (e.g. lightning, water, fire).

Improper use or transport.

Failure to follow the safety and maintenance instructions.

Other **improper treatment or modification** of the product.

After the expiration of the warranty period, you may wish to have your product repaired at your own expense.

This warranty gives you specific legal rights, and you may also have other rights which vary from state to state but this warranty is not intended to exceed statutory minimum requirements beyond what is included above.

This warranty does not limit the statutory obligation of the seller arising from any implied warranties. The period of warranty can only be extended in accordance with applicable law.

Neither the service company nor ALDI will assume any liability for data or settings stored on any returned product.

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